



Installation of Optimum+ TWAIN drivers

Run the installation files for the TWAIN application and follow the instructions to install.


Double-click on the installation file called: **“Setup_optimum+Twain_x32.exe”**.

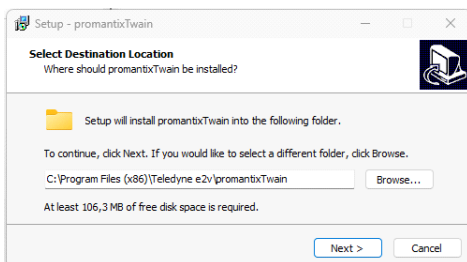


It will ask you where to install it, just leave the default installation path as it is and click on **“Next”** to start the installation.

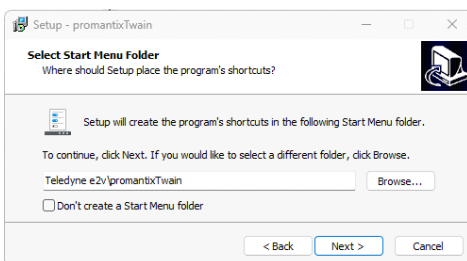


Once again just leave the default settings and click on **“Next”**.

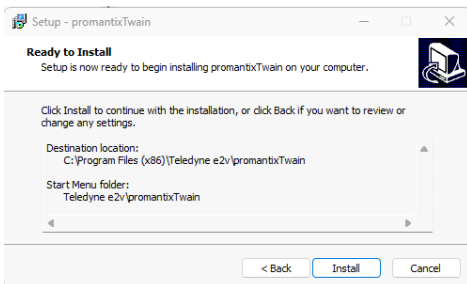


On this window only click on **“Install”** to start the installation.

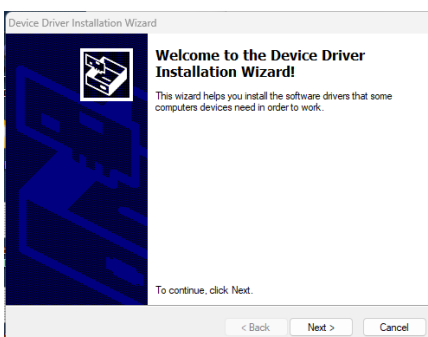


When you get this window just click on the **“Next”** buton.

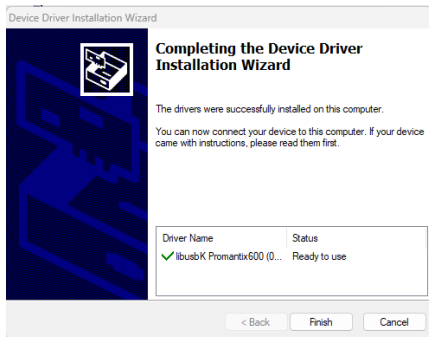


When this part is done you will see this window. Click on **“Finish”**.

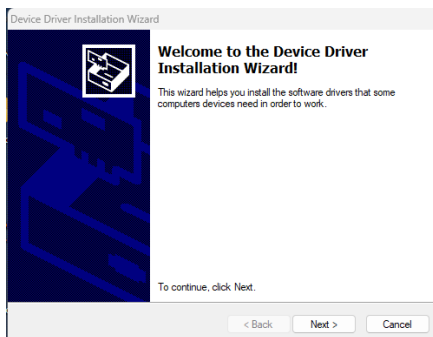

On that window that just opened, click on **“Next”**.

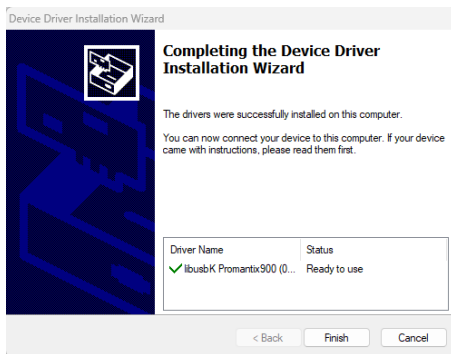

When it gets to this window, click on **“Finish”** to continue.

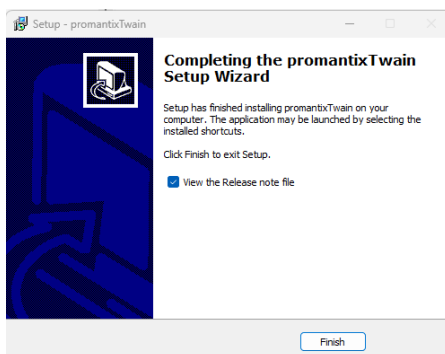

When you get this window, the installation is finish.

You can either read or not the release note. If you don't want to, just uncheck the **“View the Release note file”**.

Click on **“Finish”** to close the installation.

This must be done for EVERY computer that you want to use.

For more assistance contact **CADI** Technical Support.